





Perry L. Parke  
President/Owner


**"No Questions Asked Guarantee."**

"We unconditionally guarantee your total satisfaction with the quality and performance of all our products. If you are not completely satisfied, we will pick up the product and gladly refund your full purchase price."

ORIGINAL (WHITE)	500 PER PACK	TOTAL
	_____ @ \$4125	\$ _____
LONG (BLUE)	375 PER PACK	TOTAL
	_____ @ \$4125	\$ _____
	US Shipping	\$ <u>695</u>
	(Add sales tax if delivered in GA)	\$ _____
	<b>Total</b>	\$ _____
	Canada Shipping	\$ <u>2000</u>
	Canada Only - Free Shipping on Orders Over \$400	\$ _____
	Canadian Orders Apply GST and/ or HST	\$ _____
	(Canadian Orders ship direct from Canada) <b>Total</b>	\$ _____

# FLAPS®

Comfortable, Kind Bite Tabs

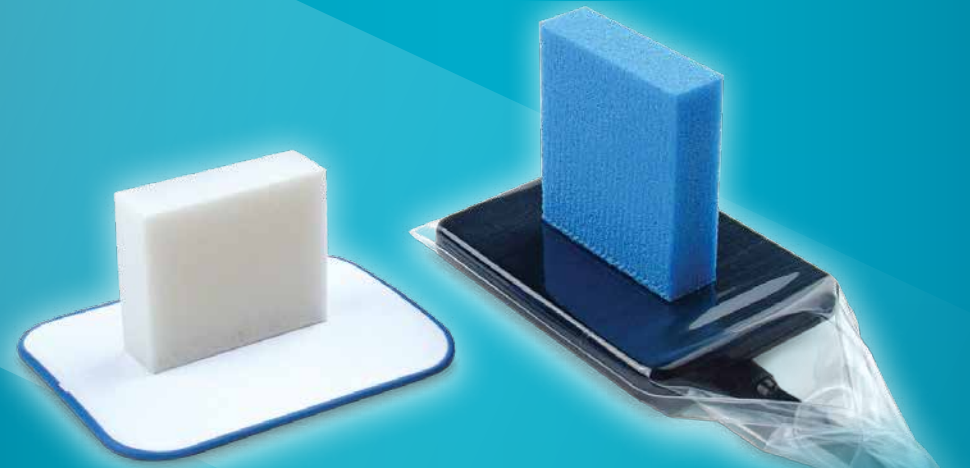

The same people who make FLAPS,  
also bring you NEODRYS, BITE-CHEK & GAZELLE.


Order Online:  
FlapsBiteTabs.com

By phone:  
800.235.1863

By fax:  
770.423.4996



**NEODRYS**  
Ultimate Parotid  
Saliva Control

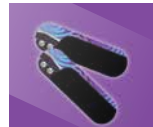

**BITE-CHEK**  
Get a Handle on  
Articulation



**GAZELLE**  
Go the Distance

“ Our patients are very happy we are switching from the cardboard stick-on-tabs. Flaps are much more comfortable!”

S. Nichols, RDH

FLAPSBITETABS.COM



Microcopy  
800.235.1863 | P.O. Box 2017 | Kennesaw, GA 30156-9017

facebook.com/microcopydental @microcopydental

Prices effective January 1, 2015

FLDT 0415



## EFFICIENT

Simple, self-sticking foam tabs attach quickly and are easy to position.

## EXACT

Cushioned tab holds film or sensor securely, allows full occlusion and provides visual reference for accurate placement and alignment.

## COMFORTABLE

Well-accepted by both children and adults. Soft cushioning, without hard surfaces, provides the ultimate in comfort and small size reduces gagging.

## DISPOSABLE

Use-and-toss convenience — extremely cost effective too!



Perfect for horizontal or vertical bitewings. Excellent for any periapical using parallel or bisecting angle technique.

Available in two sizes, Original (WHITE) and Long (BLUE).

*"These Flaps save time and are very easy to use."*



*M. Foley, DDS*

*"Flaps are very easy, simple, and fast to use. Plus, patients don't gag like they do with other devices."*

*Tamara Simac RDA*

*"I work in a pediatric dental office and the Film Flaps (Original), are great for kids. Many of my patients say they are much more comfortable."*


*Pam Melollon RDH*

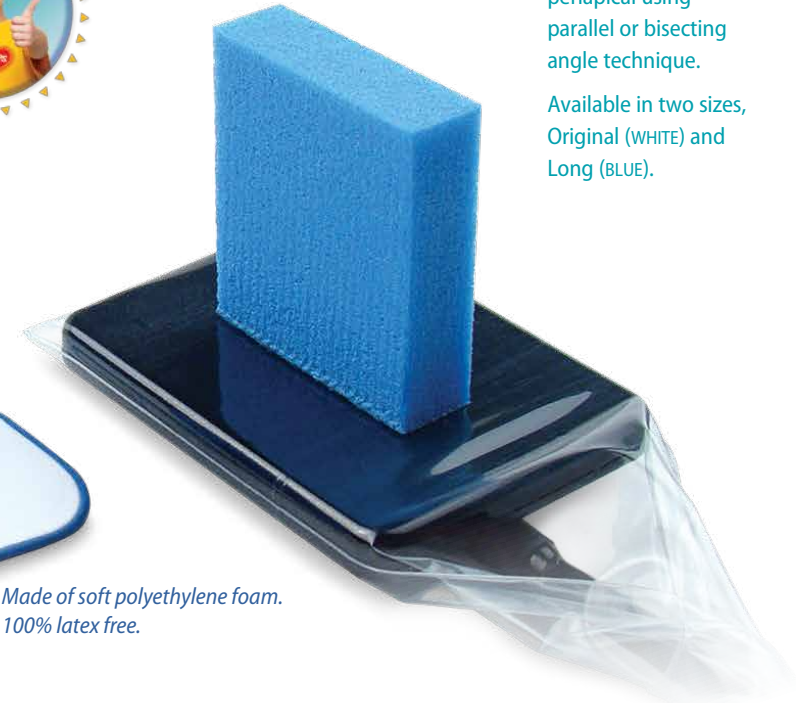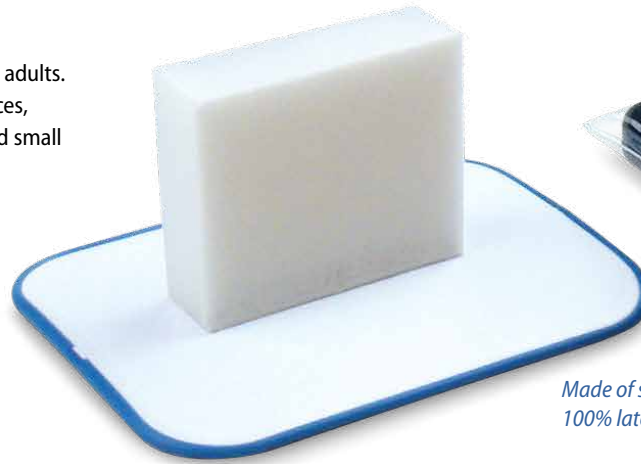


Made of soft polyethylene foam. 100% latex free.

## ORIGINAL (WHITE)


Accommodates all #0, #1, and #2 film.

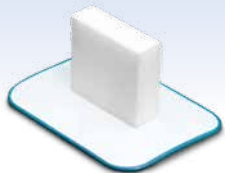


### Horizontal & Vertical Bitewings

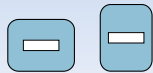


### Posterior & Anterior Periapicals

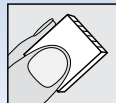

Stick the Flaps to the back of the film for the desired position.

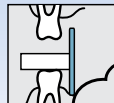


Position film. Tongue presses against back of film and teeth bite into center of the Flaps.

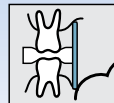


## LONG (BLUE)



Accommodates all #0, #1, and #2 sensors.

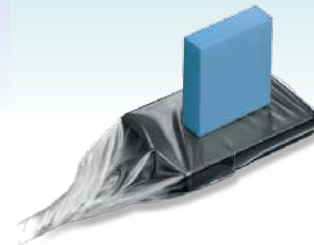

### Horizontal & Vertical Bitewings

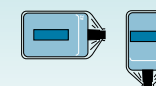


### Posterior & Anterior Periapicals

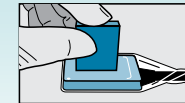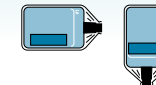

Stick the Flaps on the image side of the sensor for desired position. (See diagram)

*Required: snug-fitting protective sheaths!*

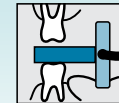

Position sensor. Tongue presses against back of sensor and teeth bite into center of Flaps. Cushion tab allows full occlusion and visual reference for image positioning.

